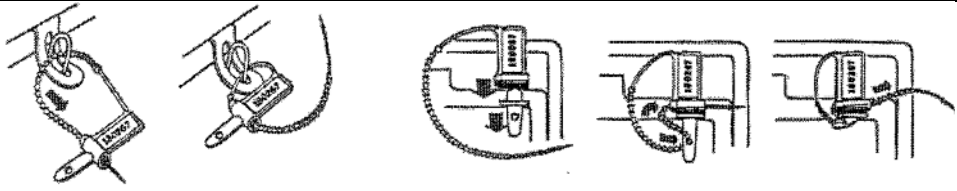


Instructions for Use: Cynch-Loks

Brand Name of Product	Cynch-Loks
Generic Name of Product	Tamper Evident Seals
Product Code Number(s)	6020 DB, 6020 GN, 6020 OR, 6020 RD, 6020 WT, 6020 YL, 6323 BK, 6323 GN, 6323 RD, 6323 RD BIO, 6323 WT, 2220 LOK
Purpose of Product	To show tampering and/or access to critical equipment and supplies.
Range of Applications for Product	Crash carts, emergency boxes, cabinets, and other areas where demonstrating tampering is critical.
Key Specifications of Product	<ul style="list-style-type: none"> Numbered or unnumbered Cynch-style allows for adjustment to application lengths Seal can be broken by hand Cannot be resealed once broken Variety of colors available 2220 LOK communicates “Clean” and “Biohazard” during transportation

Shipping & Storage	
Shipping Conditions & Requirements	N/A
Storage Conditions	N/A
Packaging Contents	100 per box
Shelf Life	N/A

Instructions for Using Product	
Description of Use(s)	Seal crash carts, emergency boxes and other areas where evidence of tampering is critical
Preparation for Use	N/A
Diagrams (drawings, pictures)	
Steps for Use of Product	<p>For Normal Apertures:</p> <ol style="list-style-type: none"> Detach bail tip from locking barrel. Pass bail through opening of item to be sealed and through tapered end of locking barrel. Pull bail tightly all the way through. Tug on seal to ensure it is securely locked. <p>For Narrow Slot Apertures:</p> <ol style="list-style-type: none"> Detach bail tip from locking barrel. Insert tab into narrow slot. Bend tab upwards & pass bail through hole of tab. Twist bail around & insert into tapered end of locking barrel. Pull bail through barrel as far as it will go.
Interpretation of Results	N/A
Contraindications of Test Results	N/A
Documentation	Follow the policy of your facility for recording date, staff and other information critical to management of critical inventory.
Special Warnings and Cautions	If a seal is broken, assume items have been tampered with or removed. Follow your facility instructions for inspection and replacement of critical inventory.
Disposal	N/A

Reprocessing Instructions	
Point of Use	N/A
Preparation for Decontamination	N/A
Disassembly Instructions	N/A
Cleaning – Manual	N/A
Cleaning – Automated	N/A
Disinfection	N/A

Drying	N/A
Maintenance, Inspection, and Testing	N/A
Reassembly Instructions	N/A
Packaging	N/A
Sterilization	N/A
Storage	N/A
Additional Information	N/A

Related Healthmark Products	N/A
Other Product Support Documents	Security Products Brochure, Security Products Price List
Reference Documents	N/A
Customer Service Contact	Healthmark Industries Company, Inc. 18600 Malyn Blvd. Fraser, MI 48026 1-586-774-7600 healthmark@hmark.com hmark.com