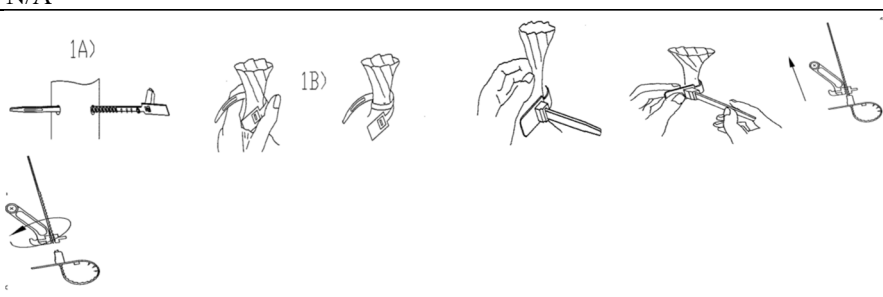


Instructions for Use: Long Cynch-Loks

Brand Name of Product	Long Cynch-Loks
Generic Name of Product	Tamper Evident Seals
Product Code Number(s)	927-15 BL, 927-15 RD, 927-07 YL, 927-11 BL, 1515 GN HLD
Intended Use	To show tampering and/or access to critical equipment and supplies.
Range of Applications for Product	Crash carts, emergency boxes, cabinets, and other areas where demonstrating tampering is critical.
Key Specifications of Product	<ul style="list-style-type: none"> • Numbered • Seal and broken by hand • Cannot be resealed once broken • Variety of colors available

Shipping & Storage	
Shipping Conditions & Requirements	N/A
Storage Conditions	N/A
Packaging Contents	N/A
Shelf Life	N/A

Instructions for Using Product	
Description of Use(s)	Seal crash carts, emergency boxes, and other areas where evidence of tampering is critical.
Preparation for Use	N/A
Diagrams (drawings, pictures)	
Steps for Use of Product	<ol style="list-style-type: none"> 1. Place strap through locking chamber with teeth towards neck of bag. 2. Pull strap manually or place hook of puller tool in slotted strap-end. Pull tightly. 3. For tight grip, place jaws of the tool over teeth and pull. The bag is now securely sealed. 4. To remove surplus strap, place jaws of tool over teeth at locking chamber and twist strap.
Interpretation of Results	N/A
Contraindications of Test Results	N/A
Documentation	Follow the policy of your facility for recording the date, staff, and other information pertinent to management of critical inventory.
Special Warnings and Cautions	If a seal is broken, assume items have been tampered with or removed. Follow your facility's instructions for inspection and replacement of critical inventory.
Disposal	N/A

Reprocessing Instructions	
Point of Use	N/A
Preparation for Decontamination	N/A
Disassembly Instructions	N/A
Cleaning – Manual	N/A
Cleaning – Automated	N/A
Disinfection	N/A
Drying	N/A

Maintenance, Inspection, and Testing	N/A
Reassembly Instructions	N/A
Packaging	N/A
Sterilization	Compatible with steam sterilization. Standard cycle: 132 °C (270 °F) for 4 minutes
Storage	N/A
Additional Information	N/A

Related Healthmark Products	N/A
Other Product Support Documents	Security Products Brochure, Security Products Price List
Reference Documents	N/A
Customer Service Contact	Healthmark Industries Company, Inc. 18600 Malyn Blvd. Fraser, MI 48026 1-586-774-7600 healthmark@hmark.com hmark.com